




**Your  
Smile is  
Priceless**

**ruby|dent**

**ENDOART®**  
FINAL TOUCH

Product Catalogue



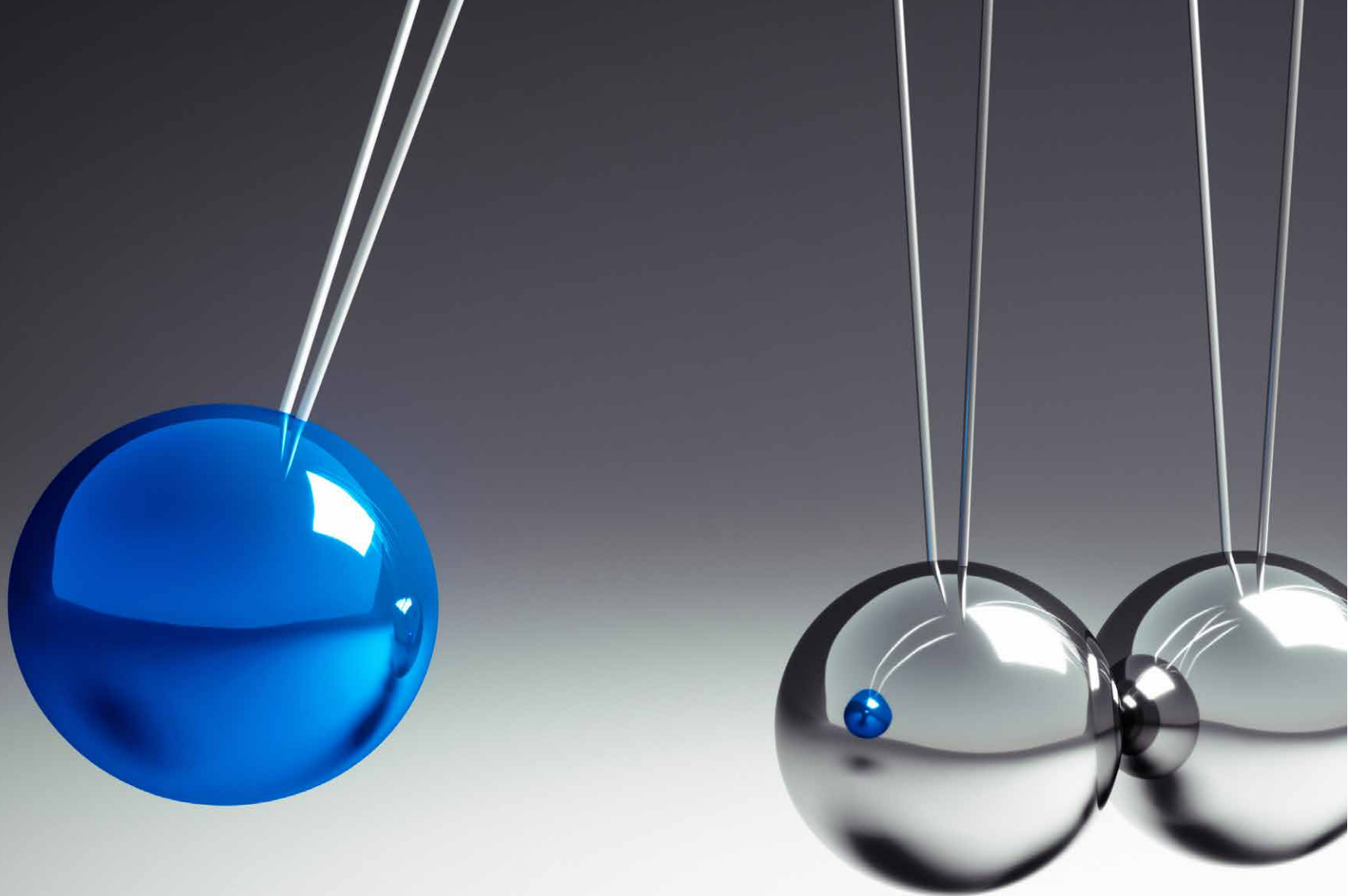

 incidental

# INDEX

## INTRODUCTION

## RESTORATIVES

RubyCap Capsule Amalgam  
RubyCap NG Capsule Amalgam  
RubyCap High Capsule Amalgam  
RubyGam Powder Amalgam  
RubyGam High Powder Amalgam  
RubyComp Nano Composite  
RubyComp P Posterior Composite  
RubyExcellence Nano Composite  
RubyFlow Flowable Composite  
RubyEtch Etching Gel  
RubyLiner  
Ruby Pit and Fissure Sealant  
RubyBond  
RubySEBond

## ENDODONTICS

RubyCal Calcium Hydroxide Paste With Barium Sulphate  
RubyCal Calcium Hydroxide Paste With Iodoform  
RubyCal LC Calcium Hydroxide  
RubyEdta EDTA Gel  
EndoArt Rotary Files  
EndoArt Reciproc Files  
EndoArt Pedo Files  
EndoArt Path File  
EndoArt K Files  
EndoArt H Files  
EndoArt Reamers  
EndoArt Barbed Broaches (Tirnerf)  
EndoArt Finger Plugger  
EndoArt Finger Spreader  
EndoArt Endomotor with Apex Locator  
EndoArt Apex Locator  
EndoArt Gutta Percha Points  
EndoArt Absorbent Paper Points  
EndoArt Fiber Post

**5** **PROPHYLAXIS**  
RubyShine Prophylaxis Paste **27**  
FlexDam Gingiva Barrier **27**

**6**  
**6** **POLISHING & FINISHING**  
**7** RubyPlaton Polishing Discs **28**  
**8** RubyPlaton Polishing Strips **28**

**8**  
**10** **WEDGES**  
**10** RubyPlaton Wooden Wedges **29**  
**11** RubyPlaton Transparent Wedges **29**

**12**  
**12** **ARTICULATING PAPERS**  
**13** RubyMark Articulating Paper Straight **30**  
**13** RubyMark Articulating Paper Horseshoe **31**

**14**  
**14** **STUDY MODELS**  
Arma Dental Study Model **32**  
Arma Dental Phantom Teeth **32**

**16**  
**16** **AUXILIARY**  
**17** RubyCrownRemover **33**

**17**  
**19**  
**20**  
**21**  
**21**  
**22**  
**22**  
**22**  
**23**  
**23**  
**23**  
**24**  
**24**  
**25**  
**25**  
**26**





### **Company**

Inci Dental was founded 1994 in Istanbul. Since then the company produce a variety of consumable materials. Since 1994, the concept of quality is established, always giving priority to customer satisfaction, the production process by making significant investments in R & D activities, cooperate with universities and Scientific and Technological Research Council of Turkey, made production as newest and unique products, in the sector so has achieved a leading place.

The company is only producer of variety materials in Turkey. Inci Dental is manufacturer of Dental amalgam filling material, dental composites filling material, endodontic hand files and rotary files, fiber posts, polishing discs, wedges, acrylics and etc. We export our products more than 45 countries in the World.

Our quality assurance approve by European Union with the certificates of ISO 9001, ISO 13485 and the EC.

### **Vision**

To be a company which is the first choice of domestic and foreign business partners.  
To produce best quality materials as only producer in the country.  
In this context, the sale of the products that we have done all, are the first brands that come to mind of the dentists.

### **Mission**

Inci dental, in dental sector;  
Who makes the difference  
Who find, research and use the right technologies in all over the world  
Who adapt quickly to technological developments  
Who produce best quality  
Who lead the innovation  
A company which holds the highest level of employee, domestic and foreign business partners satisfaction.



## RubyCap Capsule Amalgam

High copper alloy which contains a mixture of spherical and irregular particles. Predosed in self activating capsules. Gamma-2 free. It has 4 spill options as spill 1, spill 2, spill 3 and spill 5.

### Indications

Stress bearing fillings in posterior teeth (Class I and II)  
Retrograde root canal filling after apex resection

### Features

Mixture of spherical and irregular particles  
High copper alloy  
High affinity with mercury  
Firm condensation  
Smooth carving  
Excellent condensability

### Composition

45% Ag – 30,5% Sn – 24% Cu – 0,5% Zn  
45% Ag – 31% Sn – 24% Cu (Zinc Free)

### Packaging

50 capsule in a jar



## RubyCap NG Capsule Amalgam

High copper alloy which contains a mixture of spherical and irregular particles. It contains higher silver (55% Ag) than RubyCap. Predosed in self activating capsules. Gamma-2 free. It has 4 spill options as spill 1, spill 2, spill 3 and spill 5.

### Indications

Stress bearing fillings in posterior teeth (Class I and II)  
Retrograde root canal filling after apex resection

### Features

High silver / high copper alloy  
Easy to use

### Composition

55% Ag – 29% Sn – 15,8% Cu – 0,2% Zn

### Packaging

50 capsule in a jar



## RubyCap High Capsule Amalgam

High silver alloy which contains a mixture of spherical and irregular particles. Predosed in self activating capsules. Gamma-2 free. It has 4 spill options as spill 1, spill 2, spill 3 and spill 5.

### Indications

Stress bearing fillings in posterior teeth (Class I and II)  
Retrograde root canal filling after apex resection

### Features

Mixture of spherical and irregtor particles  
High silver content  
Exceptionally durable  
Easy to process  
Relatively fast hardening

### Composition

%69,2 Ag %18,6 Sn %11,9 Cu %0,3 Zn  
%69,2 Ag %18,9 Sn %11,9 Cu (Zinc Free)

### Packaging

50 capsule in a jar

### Technical Properties

A well balanced alloy for quality fillings;  
Spill 1 – 400 mg Alloy / 360 mg Mercury  
Spill 2 – 600 mg Alloy / 540 mg Mercury  
Spill 3 – 800 mg Alloy / 720 mg Mercury  
Spill 5 – 1200 mg Alloy / 1080 mg Mercury

Working Time: 4 minutes  
Carving & burnishing: 5 – 9 minutes

Conforms to:  
ISO 20749 Pre-capsulated Dental Amalgam

Compressive Strength;  
1h < 397 MPa  
24h < 437 MPa

Creep : %0,99  
Dimensional Change During Hardening :  
%0,01

RubyCap I RubyCap II RubyCap III



RubyCap NG I RubyCap NG II RubyCap NG III



RubyCap High I RubyCap High II RubyCap High III





## RubyGam Powder Amalgam

High copper alloy which contains a mixture of spherical and irregular particles. Gamma-2 free.

### Indications

Stress-bearing fillings in posterior teeth (Class I and II)  
Retrograde root canal filling after apex resection

### Features

- Mixture of spherical and irregular particles
- High copper alloy
- High affinity with mercury
- Firm condensation
- Smooth carving
- Excellent condensability

### Composition

45% Ag – 30,5% Sn – 24% Cu – 0,5% Zn

### Packaging

Boxes of 30 gr



## RubyGam High Powder Amalgam

High silver alloy which contains a mixture of spherical and irregular particles. Gamma-2 free.

### Indications

Stress bearing fillings in posterior teeth (Class I and II)  
Retrograde root canal filling after apex resection

### Features

- A mixture of spherical and irregular particles
- High silver content
- Exceptionally durable
- Easy to process
- Relatively fast hardening

### Composition

69,2% Ag – 18,6% Sn – 11,9% Cu – 0,3% Zn

### Packaging

Boxes of 30 gr




# RUBYCAP

STABILITY

HIGH STRENGTH

EXCELLENT HANDLING

NON GAMMA 2

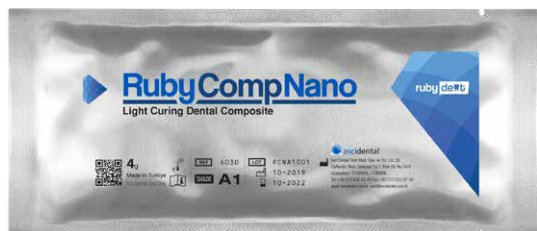

## RubyComp Nano Composite

RubyComp Nano is radiopaque nano hybrid composite system which can be the photocurable applied in the anterior and posterior field.

### Properties

- Universal Nano Composite
- It has light refraction property like a normal tooth due to the nano fillers in it
- It has high filler rate
- It can be polished high with high degree
- It has high resistance to abrasion
- It provides ease of use

### Packaging

- Refill – 4 gr in a syringe
- Kit 1 – 4 gr syringe x 6 pcs. , 5 ml RubyBond x 1 pcs. 3 ml RubyEtch x 1 pcs. and Accessories
- Kit 2 – 4 gr syringe x 8 pcs. , 5 ml RubyBond x 1 pcs. 3 ml RubyEtch x 1 pcs. and Accessories
- Shades – A1 – A2 – A3 – A3,5 – B1 – B2 – C1 – D3 – OA2 – Extra Bleach



## RubyComp P Posterior Composite

RubyComp P is micro hybrid composite system which can be the photocurable applied in the posterior field.

### Properties

- Non stick to instruments
- It has surface hardness is like real enamel
- It has high intensity
- It has a perfect consistence can be shaped/formed easily
- It has high filler rate
- It resists to chewing force
- There is not discoloration after the polymerization.

### Packaging

- Refill – 4 gr in a syringe
- Kit 1 – 4 gr syringe x 6 pcs. , 5 ml RubyBond x 1 pcs. 3 ml RubyEtch x 1 pcs. and Accessories
- Shades – A1 – A2 – A3 – A3,5



## RubyExcellence Nano Composite

RubyExcellence Nano is radiopaque nano composite system which can be the photocurable, applied in the posterior and anterior field.

### Properties

- Higher Esthetic
- Low shrinkage
- It contains Bis- GMA, UDMA, Bis- EMA resine
- Non sticky
- It has high filler rate
- It can be polished with high degree
- It has high resistance to abrasion
- It provides ease of use
- Long durability
- Non sticky

### Packaging

- Refill – 4 gr in a syringe
- Kit 1 – 4 gr syringe x 6 pcs. , 5 ml RubyBond x 1 pcs. 3 ml RubyEtch x 1 pcs. and Accessories
- Kit 2 – 4 gr syringe x 8 pcs. , 5 ml RubyBond x 1 pcs. 3 ml RubyEtch x 1 pcs. and Accessories
- Shades – A1 – A2 – A3 – A3,5 – B1 – B2 – C1 – D3 – OA2

### Technical Properties

Conforms to:  
ISO 4049 Polymer-Based Restorative Materials

- Water Solubility ( $\mu\text{g}/\text{mm}^3$ ) – 1,45  $\mu\text{g}/\text{mm}^3$
- Water Sorption <40  $\mu\text{g}/\text{mm}^3$
- Flexural Strength – 130 MPa
- Compressive Strength – 330 MPa
- Fracture Strength – 1,15 Kic
- Depth of Cure >1,5mm



## RubyFlow Flowable Composite

RubyFlow is microhybrid flowable light curing radiopaque composite.

### Properties

- Low shrinkage, low wear
- High gloss polishability
- Beautiful esthetic restorations
- Very low surface roughness
- Excellent handling
- Non sticky
- Very good marginal adaptation

### Packaging

Refill – 2 gr in a syringe  
 Kit 1 – 2 gr in a syringe x 2 pcs.  
 Kit 2 – 2 gr syringe x 4 pcs.  
 Kit 3 – 2 gr syringe x 4 pcs. , 5 ml RubyBond x 1 pcs. 3 ml  
 RubyEtch x 1 pcs. and Accessories  
 Shades – A1 – A2 – A3 – A3,5 – B1 – OA2



## RubyEtch Etching Gel

RubyEtch is an etching gel containing %37 phosphoric acid. Delivered in 3 ml. and 50 ml syringe with very thin application needles.

### Properties

- Easily Rinse
- Blue coloured
- Perfect Surface Retention

### Packaging

3 ml syringe  
 3 ml syringe x 2  
 3 ml syringe x 4  
 50 ml syringe



## RubyLiner with Glass Ionomer

Light curing glass ionomer based cavity liner to protect pulp of the tooth. It is used under direct and indirect restorations, like composites, amalgams, ceramics and another metal restorations.

### Properties

- Excellent biocompatibility
- Releases flouride
- Radiopaque
- Excellent handling characteristics

### Packaging

- 2 gr syringe
- 2 gr syringe x 2
- 2 gr syringe x 4



## Ruby Pit and Fissure Sealant

Light curing sealant with flouride releasing for occlusal surfaces. It can be used as preventive products.

### Properties

- Flouride releasing
- Radiopaque
- Easy application
- Low shirinkage

### Packaging

- 2 gr syringe
- 2 gr syringe x 2
- 2 gr syringe x 4



## RubyBond

RubyBond is light curing, single component bonding agent for the adhesive restorations. It universal and compatible with all composites.

### Properties

Low abrasion  
Elastic connection

### Packaging

5 ml in a bottle



## RubySEBond

RubySEBond is self-etching, light curing one component adhesive for the adhesive restorations.

### Properties

Very strong bonding to enamel and dentine  
No extra etching with phosphoric acid  
Long lasting restorations  
Time savings

### Packaging

5 ml in a bottle

**Your  
Smile is**  
Priceless



incidental





## RubyCal Calcium Hydroxide Paste With Barium Sulphate

Calcium hydroxide premixed paste non-setting for root canal therapy and vital pulp therapy. It comes with application tips.

### Properties

- High Radiopacity
- Easy removal
- Paste has pH > 12
- Excellent handling characteristics

### Packaging

- 2 gr syringe
- 2 gr syringe x2
- 2 gr syringe x4



## RubyCal Calcium Hydroxide Paste With Iodoform

Calcium hydroxide premixed paste non-setting for root canal therapy and vital pulp therapy. It comes with application tips.

### Properties

- Radiopacity
- Promote apexification
- Disinfects Canal
- Excellent handling characteristics

### Packaging

- 2 gr syringe
- 2 gr syringe x2
- 2 gr syringe x4



## RubyCal LC Calcium Hydroxide

Light curing calcium hydroxide for pulp capping and liner material under dental filling.

### Properties

Direct & indirect pulp capping  
Light curing  
Single component

### Packaging

2 gr syringe  
2 gr syringe x2  
2 gr syringe x4



## RubyEdta EDTA Gel

It is root canal preparation gel for lubricating and cleaning of root canal. It comes with application tips.

### Properties

Easily moving  
Antibacterial effect  
Reduce file breakage

### Packaging

3 ml syringe  
3 ml syringe x2  
3 ml syringe x4



**ENDOART®**  
FINAL | TOUCH



**Heat Treated Control Memory  
Rotary Files**



ROTARY FILES

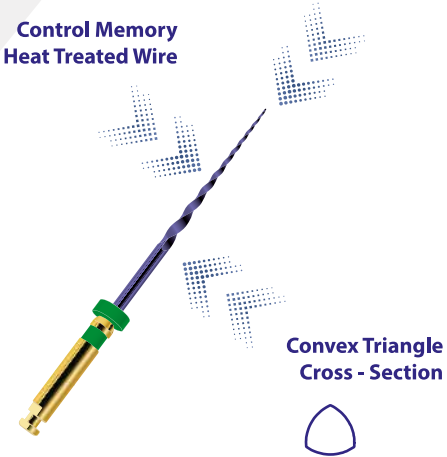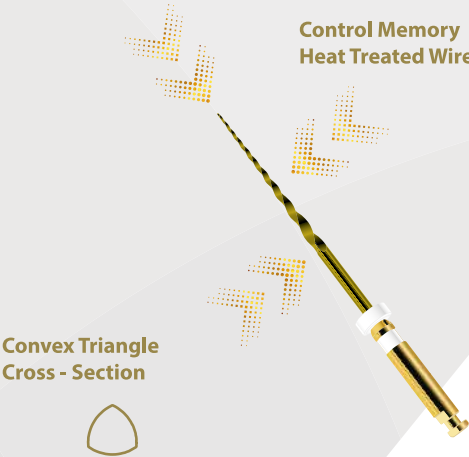

EndoArt Ni-Ti Gold/Blue Rotary Files


|                                                      |                         |                     |             |
|------------------------------------------------------|-------------------------|---------------------|-------------|
| Action Kit Sx - S1 - S2 - F1 - F2 - F3               |                         |                     |             |
| Gold                                                 |                         | Blue                |             |
| Starter File                                         | : Sx                    | Torque: 3.5 Ncm     | 300/350 Rpm |
| Shaping Files                                        | : S1 - S2               | Torque: 1.5/2.0 Ncm | 300/350 Rpm |
| Finishing Files                                      | : F1 - F2 - F3          | Torque: 2.5/3.0 Ncm | 300/350 Rpm |
| Touch Kit Sx - 15/04 - 20/04 - 20/06 - 25/06 - 30/06 |                         |                     |             |
| Gold                                                 |                         | Blue                |             |
| Starter File                                         | : Sx                    | Torque: 3.5 Ncm     | 300/350 Rpm |
| Shaping Files                                        | : 15/04 - 20/04         | Torque: 1.0/1.8 Ncm | 300/350 Rpm |
| Finishing Files                                      | : 20/06 - 25/06 - 30/06 | Torque: 2.5/3.0 Ncm | 300/350 Rpm |
| Smart Kit Sx - 16/02 - 20/04 - 25/04 - 30/04 - 35/04 |                         |                     |             |
| Gold                                                 |                         | Blue                |             |
| Starter File                                         | : Sx                    | Torque: 3.5 Ncm     | 300/350 Rpm |
| Shaping Files                                        | : 16/02 - 20/04         | Torque: 1.0/1.8 Ncm | 300/350 Rpm |
| Finishing Files                                      | : 25/04 - 30/04 - 35/04 | Torque: 2.5/3.0 Ncm | 300/350 Rpm |

| 04 Taper |                             | 06 Taper |                             |
|----------|-----------------------------|----------|-----------------------------|
| 10/04    | Torque: 1.0 Ncm 300/350 Rpm | 10/06    | Torque: 1.5 Ncm 300/350 Rpm |
| 15/04    | Torque: 1.5 Ncm 300/350 Rpm | 15/06    | Torque: 1.5 Ncm 300/350 Rpm |
| 20/04    | Torque: 1.8 Ncm 300/350 Rpm | 20/06    | Torque: 2.5 Ncm 300/350 Rpm |
| 25/04    | Torque: 1.8 Ncm 300/350 Rpm | 25/06    | Torque: 2.5 Ncm 300/350 Rpm |
| 30/04    | Torque: 2.2 Ncm 300/350 Rpm | 30/06    | Torque: 2.5 Ncm 300/350 Rpm |
| 35/04    | Torque: 2.5 Ncm 300/350 Rpm | 35/06    | Torque: 2.5 Ncm 300/350 Rpm |
| 40/04    | Torque: 2.5 Ncm 300/350 Rpm | 40/06    | Torque: 2.5 Ncm 300/350 Rpm |
| 45/04    | Torque: 2.5 Ncm 300/350 Rpm | 45/06    | Torque: 2.5 Ncm 300/350 Rpm |
| 50/04    | Torque: 2.5 Ncm 300/350 Rpm | 50/06    | Torque: 2.5 Ncm 300/350 Rpm |

Heat Treated Wire : Gold - Blue  
Length : 21mm, 25mm, 31mm  
Packaging : 6 Files





RECIPROC

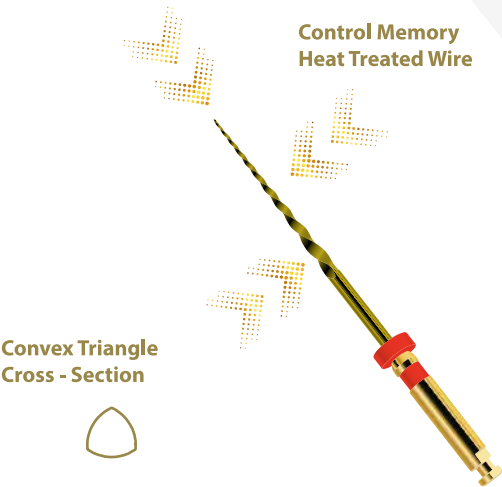

EndoArt Ni-Ti Gold  
Resiprocal Root Canal Files

Wise Reciproc Kit 20 - 25 - 35 - 45


|           |   |                  |
|-----------|---|------------------|
| Wire      | : | Gold             |
| Length    | : | 21mm, 25mm, 31mm |
| Packaging | : | 4 Files          |

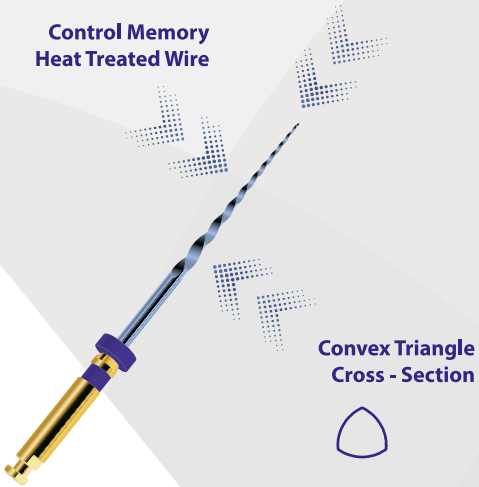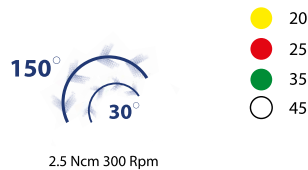

EndoArt Ni-Ti Blue  
Resiprocal Root Canal Files


Wise Reciproc Kit 20 - 25 - 35 - 45

|           |   |                  |
|-----------|---|------------------|
| Wire      | : | Blue             |
| Length    | : | 21mm, 25mm, 31mm |
| Packaging | : | 4 Files          |

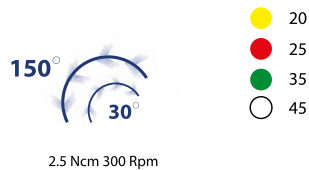




PATH FILE

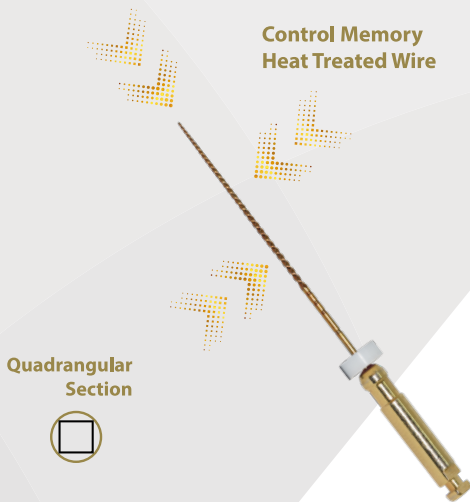

EndoArt Ni-Ti Path Rotary File

|           |                    |
|-----------|--------------------|
| Path File | 16/02              |
| Wire      | : Gold             |
| Torque    | : 1.5 Ncm 250 Rpm  |
| Length    | : 21mm, 25mm, 31mm |
| Packaging | : 6 Files          |

PEDO

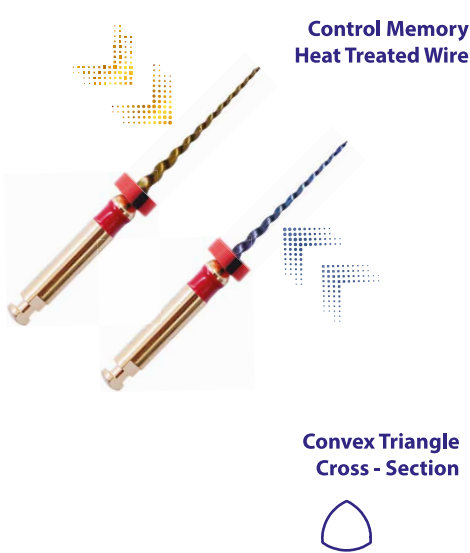

|           |                       |
|-----------|-----------------------|
| Pedo Kit  | 15/06 - 25/04 - 30/04 |
| Gold      | Blue                  |
| Wire      | : Gold - Blue         |
| Length    | : 18mm                |
| Packaging | : 3 Files             |


|       |                 |                    |
|-------|-----------------|--------------------|
| 15/06 | Torque: 1.5 Ncm | Devir: 300/350 Rpm |
| 25/04 | Torque: 1.0 Ncm | Devir: 300/350 Rpm |
| 30/04 | Torque: 2.0 Ncm | Devir: 300/350 Rpm |

## EndoArt Endomotor with Apex Locator

Cordless endomotor with apex locator function.

### Properties

- 10 Memory Settings
- 3 Rotation Mode; Clockwise/Anti Clockwise and Reciprocating
- Max 800 Rpm and 4.0 N.Cm
- Micro USB Charging
- 16:1 Contra Angle
- Auto Start in the root canal
- Auto Reverse
- 4 working mode
- Large OLED Display

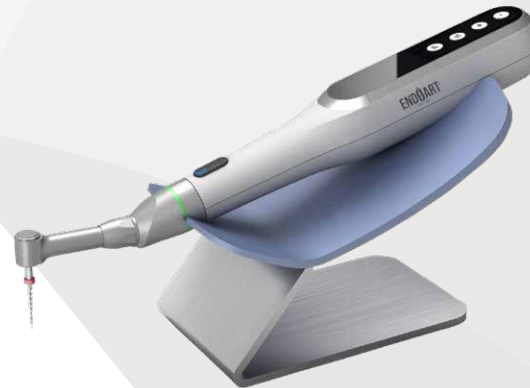

**ENDOART**  
FINAL TOUCH

## EndoArt Apex Locator

High Accuracy Apex Locator

### Properties

- Advanced Multi Frequency Tech
- Adjustable Reference Position
- Auto Calibration





## EndoArt Gutta Percha Points

Gutta Percha root canal filling material for easy fillings.  
It has variable taper and sizes.

### Properties

02/04/06 Taper Options  
Excellent Radiopacity  
Length Marked  
Hand Rolled




## EndoArt Absorbent Paper Points

Absorbent Paper Point to dry root canals.

### Properties

Highly Absorbent  
Sterile Packaging  
Variety Sizes and Tapers  
Length Marked  
Hand Rolled

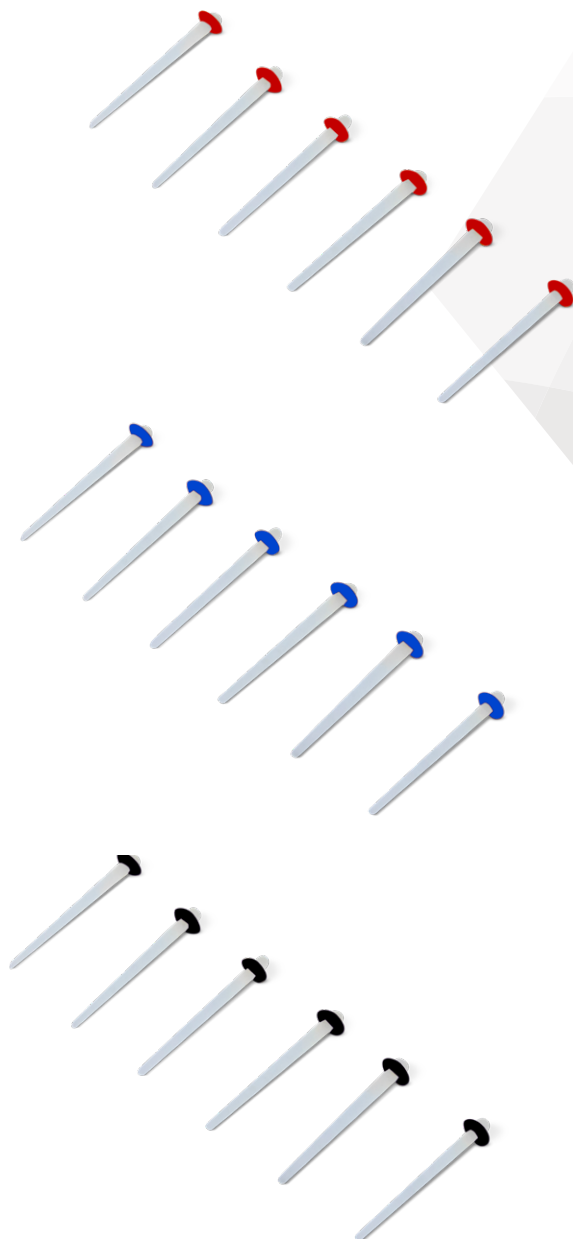

## EndoArt Fiber Post

Universal resin glass fiber post system for esthetic and metal free restorations. It has different type of diameters which are 1,2 mm, 1,4 mm, 1,6 mm and 1,8 mm.

### Properties

- Radiopaque
- 4 drills available inside of the kit
- Conical tip
- High resistance
- Colour coded
- It is compatible with adhesive resin cements

### Packaging

- EndoArt Fiber Post Kit – 5 pcs. 1,2 mm posts, 5 pcs. 1,4 mm posts, 5 pcs. 1,6 mm posts, 5 pcs. 1,8 mm posts, 4 drills
- EndoArt Fiber Post Refill – 6 pcs. 1,2 mm posts
- EndoArt Fiber Post Refill – 6 pcs. 1,4 mm posts
- EndoArt Fiber Post Refill – 6 pcs. 1,6 mm posts
- EndoArt Fiber Post Refill – 6 pcs. 1,8 mm posts



## RubyShine Prophylaxis Paste

RubyShine Prophylaxis Paste is specially formulated of polishing and cleaning agents which has nano particuls in it. Available in mint, strawberry, bubblegum and green apple flavor.

### Packaging

30 gr in a bottle  
150 gr in a bottle  
340 gr in a bottle



## FlexDam Gingiva Barrier

Light curing resin-based gingiva barrier to protect soft tissue.

### Properties

Easily removable after treatment  
It minimizes heat and soft tissue sensivity during curing  
Blue shade  
5 second curing

### Packaging

2 gr syringe  
2 gr syringe x2  
2 gr syringe x4

## RubyPlaton Polishing Discs

Has different kind of discs which has different grins and different sizes. It allows excellent results in polishing. Flexible and easy to use

### Properties

10 mm blue, yellow, white, green, brown  
12,7 mm blue, yellow, white, green, brown  
14 mm blue, yellow, white, green, brown

- **White** – Extra Fine for Polishing
- **Yellow** – Fine for Finishing
- **Green** – Medium for Countoring
- **Blue** – Coarse for Gross Reduction
- **Brown** – Extra Coarse for Gross Reduction

### Packaging

80 pcs discs and 1 mandren in a plastic box  
240 pcs discs and 4 mandren in a plastic box



## RubyPlaton Polishing Strips

Has different kind of strips which has different grin sizes. They are produced combining two different types of grins.

### Properties

It allows excellent results in polishing. Flexible and easy to use

- **White** – Extra Fine for Polishing
- **Yellow** – Fine for Finishing
- **Green** – Medium for Countoring
- **Blue** – Coarse for Gross Reduction

### Packaging

75 pcs strips in a plastic box  
80 pcs discs and 1 mandren + 75 pcs strips in a plastic box  
240 pcs discs and 4 mandren + 75 pcs strips in a plastic box





## RubyPlaton Wooden Wedges

RubyPlaton wooden wedge it has anatomic form for you to make excellent filling.

### Properties

- Superthin short
- Thin short
- Thin long
- Medium short
- Medium long

### Packaging

100 pcs. wedges in a plastic box  
200 pcs. wedges in a plastic box  
500 pcs. wedges in a plastic box



## RubyPlaton Transparent Wedges

RubyPlaton transparent wedge is formed by special transparent polymer for excellent filling.

### Properties

Not easy to break long lived with high quality and durability  
Easy to use  
Anatomic form  
Small size - Medium size

### Packaging

100 pcs in a plastic box



## RubyMark Articulating Paper Straight

Great sensitive articulating paper to get simple marking. It is straight form on the 40µ, 80µ and 200µ.

### Properties

- Non sticky
- Smooth surface
- Double sided options.

### Packaging

- 40µ – Blue/Blue Straight 144 pcs. paper
- 40µ – Red/Red Straight 144 pcs. paper
- 40µ – Blue/Red Straight 144 pcs. paper
- 80µ – Blue/Blue Straight 144 pcs. paper
- 80µ – Red/Red Straight 144 pcs. paper
- 80µ – Blue/Red Straight 72 pcs. paper
- 80µ – Blue/Red Straight 144 pcs. paper
- 200µ – Blue/Blue Straight 300 pcs. paper
- 200µ – Red/Red Straight 300 pcs. paper



## RubyMark Articulating Paper Horseshoe

Great sensitive articulating paper to get simple marking. It is horseshoe form on the 80µ.

### Properties

- Non sticky
- Smooth surface
- Double sided options.

### Packaging

80µ – Blue/Red Horseshoe 72 pcs. paper



### Arma Dental Study Model


Dental study model for students. It shows great and real-like experience during application.

#### Properties

- Hard Plastic jaw base
- Non-breakable
- The students thus practice under “human” conditions
- Custom made screws for tooth fixing
- Natural anatomy and occlusion

#### Packaging

1 piece dental study model with screwed teeth

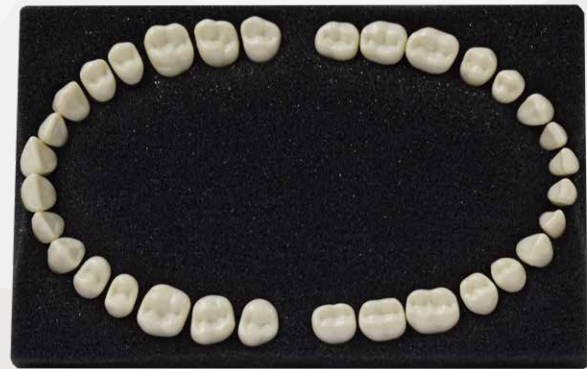

### Arma Dental Phantom Teeth

Dental study model teeth for students. It shows great and real-like experience during application. Teeth are made from melamine material.

#### Properties

- Custom made screws for tooth fixing
- Natural anatomy and occlusion

#### Packaging

1 set dental study model teeth and screws with teeth chart



## RubyCrownRemover

It is practical crown remover for easy procedure. It removes crown using the way of impact method with arc.

### Properties

Can be sterilized at 134 °C

3 types of remover tips

Different types of impact





**İnci Dental Tıbbi Malz. San. ve Tic. Ltd. Şti.**

Defterdar Mah. Otakcılar Cad.

E Blok Sit. (Has On Ticaret Merkezi) No:74 Eyupsultan / İstanbul / Turkey

**Tel : +90 212 632 63 45 Fax : +90 212 632 97 56**

[www.incidental.com.tr](http://www.incidental.com.tr) [international@incidental.com.tr](mailto:international@incidental.com.tr)